



BRODA[®]
Comfort is Essential

SYNTHESIS
TRANSPORT CHAIR

SYNTHESIS TRANSPORT CHAIR V4 WC19

COMFORT AND SAFETY ON THE GO

Providing the comfort and functionality BRODA is known for, the Synthesis V4 mobile transport chair supports patients through multiple levels of care.



- 1 INFINITELY ADJUSTABLE SEAT TILT
- 2 INFINITELY ADJUSTABLE BACK RECLINE (TRENDELENBURG CAPABLE)
- 3 HEIGHT ADJUSTABLE ARMRESTS
- 4 SWING AWAY REMOVABLE ARMS
- 5 ADJUSTABLE WINGS FOR UPPER LATERAL SUPPORT
- 6 INDEPENDENT ELEVATING, ARTICULATING, LENGTH ADJUSTABLE LEGREST
- 7 SPLIT FLIP DOWN FOOTREST
- 8 WC-19 FACTORY INSTALLED COMPONENTS (SEE SPECIFICATIONS)

Shown with standard rear casters and upgraded twin wheel front casters.

PATIENT BENEFITS



Lower center of gravity provides industry leading safety and stability. (Do not use full tilt during transport)



Lay flat feature facilitates supine transfers, reducing the risk associated with falls and injuries. (Do not recline during transport)



Standard 6" rear multi-directional locking casters provide ease of maneuverability.



Reduces the cost of non-emergent transport. Front and rear securements eliminate the need for additional transport equipment.

CAREGIVER AND FACILITY BENEFITS

- Optimum safety and functionality withstands severe agitation and involuntary movements.
- Tilt and recline provides ease of use for frequent re-positioning, helping improve patient outcomes.
- Enhanced positioning and innovative features allow for fewer transfers, promoting safe patient handling.
- Proper seating alignment decreases patient fatigue and postural deviations reducing the use of restraints.
- Quality construction provides superior durability reducing equipment maintenance and replacement costs.

BRODA[®]
Comfort is Essential

SYNTHESIS TRANSPORT CHAIR V4 WC19

COMFORT TENSION SEATING® PROVIDES

- Long term seating comfort.
- Conforming support for proper pelvic positioning.
- Improved air flow, decreasing heat and moisture build up.
- Pressure redistribution that aids in maintaining skin integrity.
- Accommodation for conditions such as kyphosis and lordosis.
- Decreased pressure, improved blood flow and oxygenation of the tissue.



CTS Strapping flexes for individualized fit & superior comfort.

SPECIFICATIONS

Frame	16 Gauge tubular steel frame
Seat & Backrest	BRODA Comfort Tension Seating®
Degrees of Tilt	40° (13° DURING TRANSPORT)
Degrees of Recline	90° (Trendelenburg capable)
Standard Seat Width	16", 18", 20" (40cm, 45cm, 50cm)
Custom Width	22", 24" (55cm, 60cm)
Standard Seat Depth	17", 18.5" (43cm, 46cm)
Standard Seat Height	20.5" (51cm)
Front Casters	5" Front total lock (12.5cm)
Rear Casters	6" Multi-directional locking (15cm)
Mag Wheel Upgrade	14" or 20" Rear mag or, 14" mid-mount wheels (35cm or 50cm)
Front Caster Upgrade	5" Aluminum front casters(rear mags only)
Weight Capacity	350lbs / 159kg
Standard Padding (In Accent Color)	Headrest cushion, shoulder bolsters, lower side pads, upper side pads, sole pads
Standard Padding (In Black Vinyl)	Seat and back pad with cover, calf pad
WC-19 Factory Installed Components	Front and rear securements, seatbelt, footrest tether, option of front crash stabilizer with standard front casters, or aluminum or twin wheel front caster

UPHOLSTERY OPTIONS




* In honor of our veterans service, camo options now available.

COMMON ACCESSORIES

- Aluminum front casters
- Twin wheel front casters
- Huntington's padding package (HSP)
- Padded armrests
- Full wrap vinyl seat and back covers
- BRODA app (additional positioning package)

FRAME OPTIONS

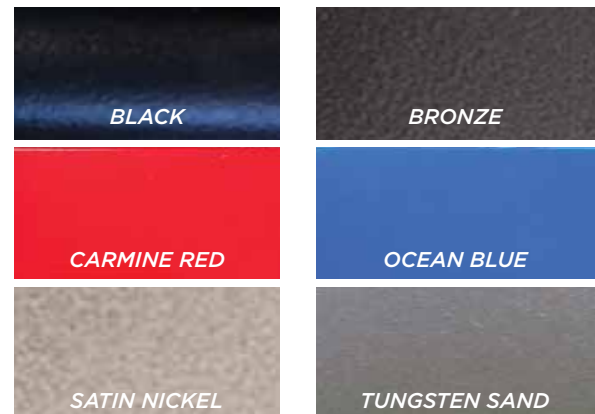



HEALTHCARE DIVISION
info@brodaseating.com
p. 1-844-552-7632

Brodaseating.com
©Broda USA, INC.
7/17